




MagCore[®] Viral Nucleic Acid Extraction Kit

For extracting Viral DNA/RNA from serum, plasma, cell-free body fluids by using the MagCore[®] HF16 System of RBC Bioscience Corp.

Cat. No. MVN400-01 // MVN400-02

Cartridge Code: 201

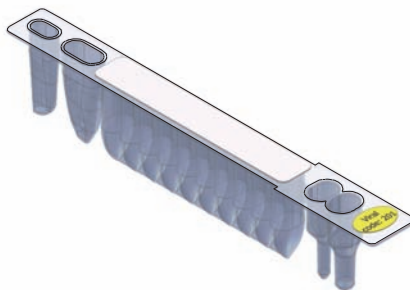

MagCore HF16 Automated DNA/RNA Purification System



MagCore® HF16 System

Reliable Nucleic Acid Purification

MagCore HF16 System is a simple, fast and cost-effective technique to automatically purify nucleic acids from a diverse range of sample sources. With the pre-programmed protocols and magnetic bead based reagent cartridges the system provides consistent and stable nucleic acid purification for every busy laboratory.

Advantage

Touch Button operation means simplicity

Rapid Results in about 30 minutes for DNA purification


Flexible Design means single samples can be run with no waste

Batch Purify up to 16 samples at one time

Diverse Purification kits and protocols available


Economic Pricing of reagents and instrument

MagCore HF16 Introduction

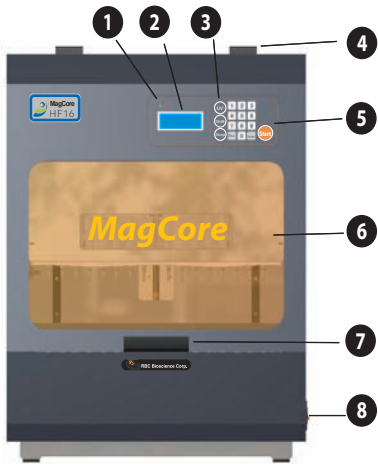


- 1. Indicator LCD
- 2. LCM Display Screen
- 3. UV Key
- 4. Aluminum Hinge
- 5. Keypress Panel
- 6. Protective Acryl Board
- 7. Door Pull
- 8. Power Switch
- 9. UV Lamp
- 10. Door-Open Detector
- 11. Door Pry
- 12. "MagCore" LCD Window
- 13. Syringe Unit
- 14. Operation Area

Instrument Supplied Accessories

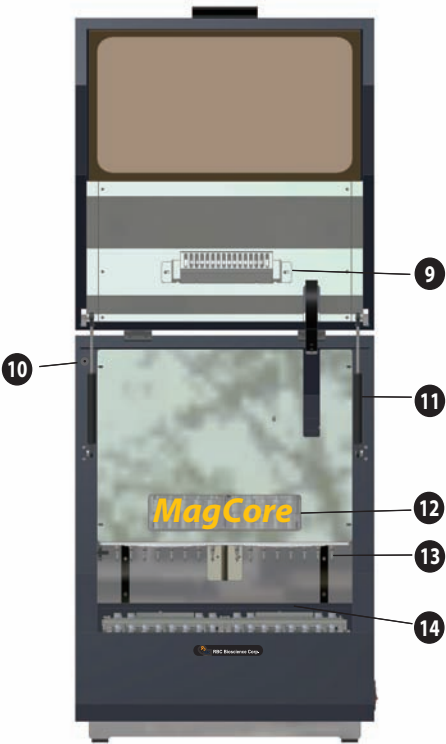

Cartridge Rack X 2

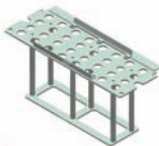

T-Rack X 2



Empty Cartridge X 20



Tip Holder X 50



Tip X 50





Sample/Elution Tube X 50

MagCore HF16 Specification

| | |
|---------------------|--|
| Model | HF16 |
| Processing Capacity | 16 samples per batch |
| Protocols V1.01 | Programs built-in. Whole Blood (200/400ul), Viral Nucleic Acid (DNA/RNA), Tissue Genomic DNA, Plant Genomic DNA |
| System Method | Cellulose coated magnetic beads |
| Sample Volume | 200/400/1,200µl (2,400µl available on request with optional 2.4ml Rack) |
| Elution Volume | 60/100/150/200 µl |
| Yield | Average 6µg Genomic DNA / 200µl human whole blood |
| Purity | O.D. 1.8-2.0 $A_{260}/_{280}$ ratio |
| Processing Time | 30-70 minutes (depends on sample type and method) |
| System Components | 1. Pipetting Unit: Dispensing, transferring, X-Z two axis movements. 2. Electric Control: Internal microprocessor 3. UV Light: Power 5W, life duration >1,000 Hrs 4. Heating Block: RT~100°C 5. Display Screen: 3 inches LCM Screen with keypress Panel 6. Accessories: T- Rack, Cartridge Rack , Empty Cartridges, Tip/Holders , Microcentrifuge Tube, Syringe O-ring and Grease |
| Power Supply | Voltage: AC 100-240V; Frequency: 50/60Hz; Power Consumption >1.0KVA |
| Dimension | W550 X D660 X H680 (mm) / W21.65 X D25.98 X H26.77 (inches) |
| Net Weight | 70kg / 155lbs |

LCM Display Screen and Keypress Panel

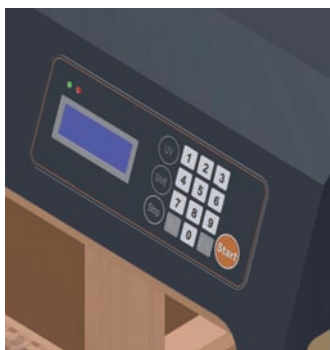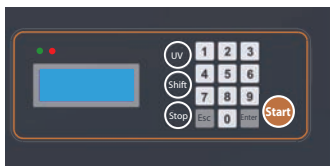

MagCore HF16 provides a touch button keypad for you to communicate with your instrument. The number keys give users to input Cartridge Code which specified on the Reagent Cartridge to run built-in program automatically, or to select desired volumes at start-up sample and final elution steps.

LCM Display Screen which with sapphire background, shows clear information that the MagCore HF16 is currently performing.

Special function of buttons described as:


- UV** UV button is designed to turn on UV Lamp while your MagCore HF16 needs to be sterilized after running contagious samples. Options with 30 / 60 mins.
- Shift** Shift button is designed for calibrating position accuracy at MagCore HF16 System, **DO NOT** execute this program without trained engineer aside.
- Stop** Stop button is designed for going back to "Stand-By" menu at program completed, or for **EMERGENCY STOP** use.

Operation Procedures



1


Turn on power switch located at the right hand side of MagCore HF16. A beep sound can be heard and the machine is on "Stand-By" condition.

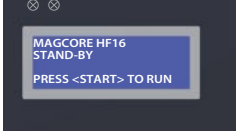

2


LCM display screen shows "Stand-By" as left picture.

Note:
At this step, press **UV** button, to perform UV decontamination

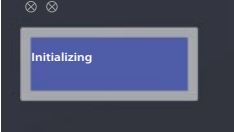


3

After press **Start** button, machine runs program of calibration, initialize to move all axis to original factory positions.

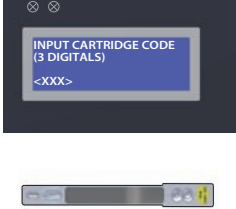


4

Input Cartridge Code : **201** to run Viral Nucleic Acid (DNA/RNA) protocol. Cartridge Code is shown on your Reagent Cartridge and the cover of user manual.

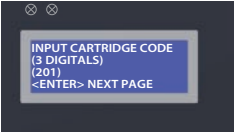

5



Confirm your input code again and press **Enter** to next page for sample volume selection.

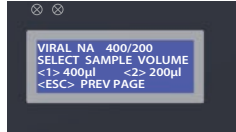

6

Viral Nucleic Acid 400/200 protocols installed.


Press **1** for 400ul sample
Press **2** for 200ul sample
Press **Esc** to re-enter cartridge code

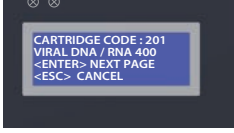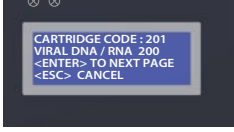7

Press **1** for 400ul sample
Press **Enter** to next page
Press **Esc** to **cancel** whole procedures and back to "Stand-By" menu.

Press **2** for 200ul sample
Press **Enter** to next page
Press **Esc** to **cancel** whole procedures and back to "Stand-By" menu.











8

At this step prepare racks to operation area. (see Step 9-18)


After racks loaded then press **Enter** to next page for elution volume selection.

Operation Procedures

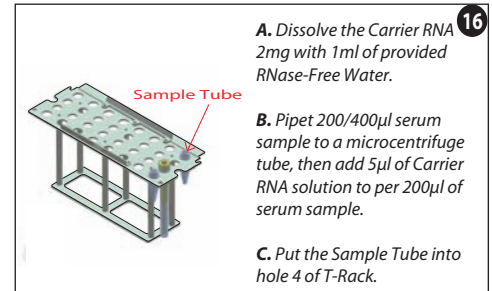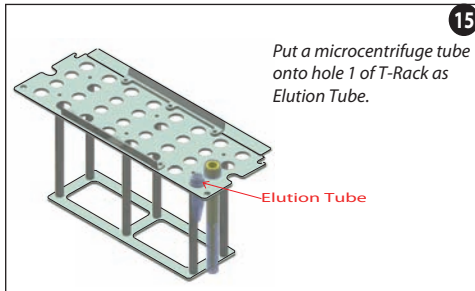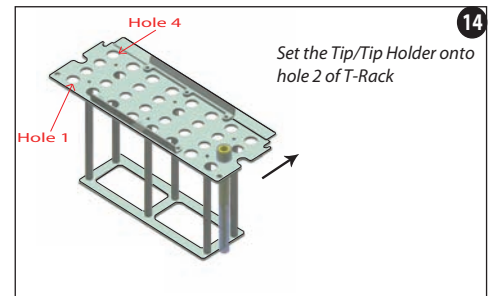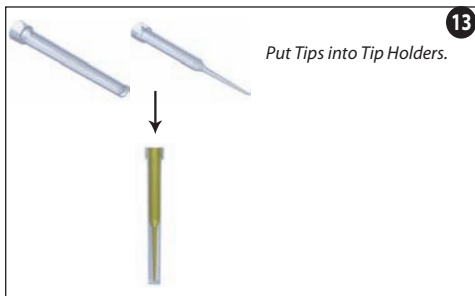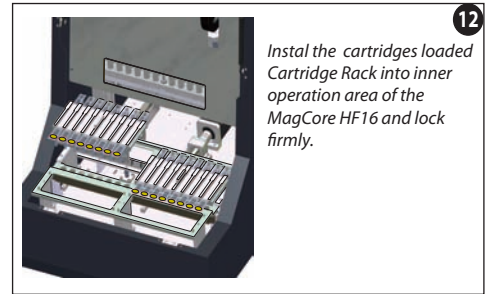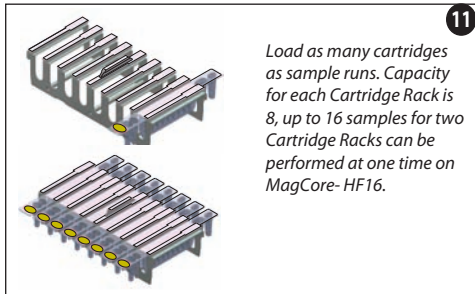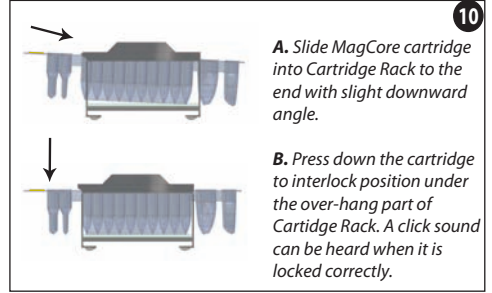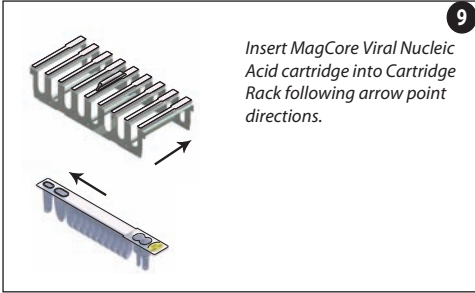


Note:
Store the Carrier RNA at 4°C for long term storage.

Operation Procedures

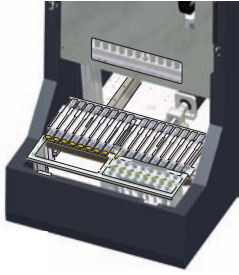


17

Install the Tip/Holder and Sample/Elution Tubes loaded T-Rack into outer operation area of the MagCore HF16 and lock firmly.



18

Pull down the front door to interlock position and then press **Enter** to next page for elution volume selection.

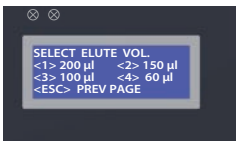

19


Select final elution volume

<1> 200µl

<2> 150µl


<3> 100µl

<4> 60µl

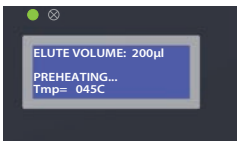


20

MagCore HF16 in process of selected protocol at this step. The Green Indicat LCD lights up and Heating Block starts to heat up to 60°C for Lysis Step.



21


During MagCore HF16 is under program running, the "MagCore" LCD lights up at all times. **DO NOT** open the door at this moment, it causes emergent stop and you might lose your samples by machine interruption.

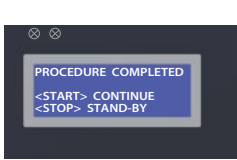

22

While program finished, a beep sound can be heard. and green Indicate LCD light went out.


LCM displays as left picutre.



23

Lift up the front door to the uppermost position, and take out Elution Tube with purified viral DNA/RNA, direct for downstram applications in PCR and RT-PCR, or store at -20°C for long term stroage.

Discard all disposables under proper procedures.

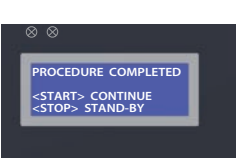

24


After clean up disposables. you can choose to

Start: run the same protocol again

Stop: to "Stand-By" menu for UV sterilization or power shut down.


UV Sterilization

1




Turn on power switch located at the right hand side of MagCore HF16. A beep sound can be heard and the machine is on "Stand-By" condition.

2

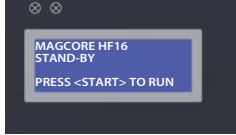


LCM display screen shows "Stand-By" menu as left picture.

3

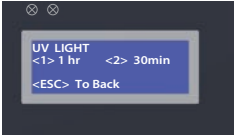

Press **UV** button to proceed UV sterilization with 30/60 mins of time options.
Esc: to "Stand-By" menu


4



UV lamp activates and count down with selected time at minutes.
MagCore HF16 will go back to "Stand-By" menu once UV sterilization procedure completed.


Door-Open Interruption / Emergent Stop

1




If the front door was opened during program performing, MagCore HF16 stops all mechanical movements immediately due to safety consideration.

2

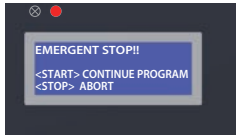


LCM display screen shows warning sign, Red Indicate LCD lights up and a beep sound can be heard.
Pull down the front door, and select
Start : continue program.
Stop : abort program and back to Stand-By" menu.

3



In case of emergency, press "STOP" button at anytime. LCM display screen shows warning sign, Red Indicate LCD lights up and a beep sound can be heard.

4



After machnical dysfunction eliminated, select
Start : continue program.
Stop : abort program and back to "Stand-By" menu.



Kit Contents

| MVN400-01 | |
|-------------------------------------|----|
| <i>Pre-filled Cartridge Reagent</i> | 36 |
| <i>Pipette Tip</i> | 36 |
| <i>Tip Holder</i> | 36 |
| <i>Microcentrifuge Tube</i> | 72 |
| <i>Tube Cap</i> | 36 |
| <i>Carrier RNA (2mg)</i> | 1 |
| <i>RNase-Free Water (10ml)</i> | 1 |

| MVN400-02 | |
|-------------------------------------|-----|
| <i>Pre-filled Cartridge Reagent</i> | 96 |
| <i>Pipette Tip</i> | 100 |
| <i>Tip Holder</i> | 100 |
| <i>Microcentrifuge Tube</i> | 200 |
| <i>Tube Cap</i> | 100 |
| <i>Carrier RNA (2mg)</i> | 2 |
| <i>RNase-Free Water (10ml)</i> | 1 |



RBC
BIOSCIENCE

RBC Bioscience Corp.

RBC Bioscience is a leading R&D based life science solutions provider.
For more information on our extensive and innovative life science range please visit our web site.